



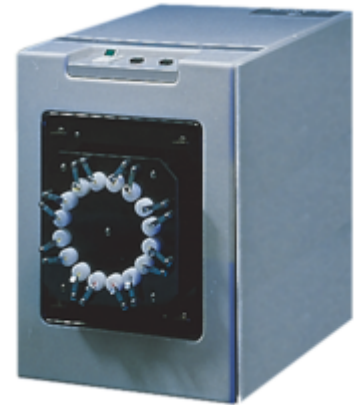

Kryo Series precision controlled-rate freezers

Kryo 360 - 1.7/3.3 - fully featured freezers for critical cell lines

- High reliability features
- Horizontal or vertical freezing
- Validation - consistent results
- Extended hold time - safe
- Multiple protocols - flexible
- Compact - less lab space
- Menu driven - simple to use

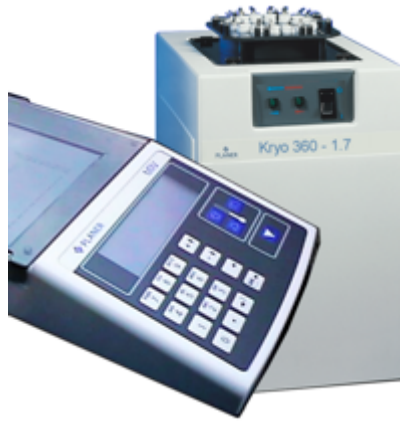

The **Kryo 360** incorporates all of the critical features expected from a high class biological freezer. The $-180\text{ }^{\circ}\text{C}$ end temperature **ensures sample integrity during transfer to storage** whilst the flexibility of the system is ideal for the more complex and demanding protocols associated with the most advanced cryopreservation techniques.

The high capacity LNP4 active nitrogen pump offers both faster cooling rates and a large reservoir, offering the **reassurance of an extended hold time** at the protocol end temperature. The system sample capacity is sufficient for the busiest laboratory and the state of the art compact design will enhance the most modern facility.

The MRV controller system has been created to offer **multiple protocols** whilst remaining simple to programme and operate and both during and after a run offers the widest range of displayed information, alphanumerically and graphically via the easy view display and as a print out on the integral full view printer.

Validation is a high priority and the MRV offers password controlled access on multiple user levels, time and date stamping, programme preview and verification before running and data storage for the last 10 runs for subsequent printing. User calibration with associated hard copy is featured and PC connection compatible with the comprehensive **Delta T™ software** application is standard.

The system has been fitted with **numerous safety features**. These help protect against power failure and PC failure when running with software. Processor or system problems are controlled and the system restarts automatically to protect samples.

For example, all control and data systems are separated, the controller can be removed from the operating freezer with no loss of programme integrity; data storage and run processing are completed on independent isolated electronic systems.

Kryo 360 Features

- Controller displays demand, sample and chamber temperatures, programme stage and current temperature graphic
- Menu driven controller, simple to programme and operate
- Horizontal or vertical operation
- Compact design
- Standard operating features:-
 - Start above ambient
 - Controlled heating
 - Data Printing (integral printer)
 - Comms port for PC connection
 - Fast cooling rates
 - Multiple safety features

SPECIFICATION OVERVIEW

- Chamber volume: 1.7 or 3.3 litres
- Straw capacity: 60 x 0.25 ml or 45 x 0.5ml in 15 positions radially located
- Ampoule capacity: (30 or 60) x 2 ml in 15 positions radially located
- Lower temperature limit: $-180\text{ }^{\circ}\text{C}$
- Cooling rates: -0.01 to $-50\text{ }^{\circ}\text{C}/\text{min}$
- Controlled heating rates: 0.01 to $10\text{ }^{\circ}\text{C}/\text{min}$
- System controller: MRV
- Delta T™ Lite Software: Upgradable to Delta T™ or Delta T™-iQ

SYSTEM REQUIREMENTS

5 psi or 22 psi Liquid Nitrogen Supply

Kryo Series precision controlled-rate freezers

Kryo 360/1.7/3.3 - fully featured freezers for critical cell lines

System Specification:

• Range	+40 °C to -180 °C
• Heating rate	0.01 °C/min to 10 °C/min
• Cooling rate	-0.01 °C/min to -50 °C/min
• Controller accuracy	±0.3 °C measured on a hold at 0 °C
• Storage temperature	-10 °C to +50 °C
• Storage humidity	5 % to 95 % relative humidity non-condensing
• Operating temperature	5 °C to 40 °C
• Operating humidity	5 % to 90 % relative humidity non-condensing

Controller Specification:

• Dimensions	80 mm high x 220 mm wide x 350 mm deep
• Weight (approx.)	2.6 kg
• Display	240 x 64 LCD with CCFL backlight
• Printer	320/640 dot thermal printer
• Keypad	20 key membrane keypad
• Programmable Cooling Rate Range	-0.01 °C/min to -99.9 °C/min
• Number of profiles	10
• Steps per profile	32
• Number of stored runs	10

Chamber Specification:

• Weight (kg)	14.4(1.7 L chamber) or 14.7(3.3 L chamber)
• Chamber dimensions	1.7 L chamber: Internal 200 mm x 150 mm diameter External 450 mm high x 300 mm wide x 420 mm deep 3.3 L chamber: Internal 400 mm x 150 mm diameter External 450 mm high x 300 mm wide x 420 mm deep
• 0.25 ml straws	60 (1.7 L and 3.3 L chamber)
• 0.5 ml straws	45 (1.7 L and 3.3 L chamber)
• 2 ml ampoules	30 (1.7 L chamber) or 60 (3.3 L chamber)
• 50 cc blood bags	-
• 250 cc blood bags	-
• 500 cc blood bags	-
• Power requirements (including MRV Controller)	115V~ 50/60Hz 600VA / 230V~ 50/60Hz 600VA

Recommendation of additional equipment:

• 5 psi System	System Pump - LNP4-C System Dewar - MVLAB 30
• 22 psi System	System Cylinder - MVEUROCYL230SB Phase Separator - MVPHASE
• 22 psi System (alternative)	Vacuum Jacketed Pipe Work System Phase Separator - MVPHASE

Planer plc, Windmill Road, Sunbury, Middlesex, TW16 7HD, United Kingdom.

Tel: +44(0)1932 755000 Fax: +44(0)1932 755001 email: sales@planer.co.uk website: www.planer.com

Planer reserves the right to change specifications without notice. All third party trademarks acknowledged. © 2010 Planer plc