



Freezer control with batch option

DeltaT Software

Deploy one profile to several freezers with our DeltaT v7 software.


Each freezing run produces a detailed, date and time stamped summary, which can automatically be stored to a secure location over a local area network.

Password protection and tamper proof files helps compliance with 21 CFR Part 11.

Profiles can be created with up to 32 different stages (heat / cool / hold) combined with real time updates of both sample and chamber sensor readings.

Email notifications during the freezing run, as well as the audible alarm of the freezer, alert users of any significant deviations from the selected profile.

- Email alerts and real time profile updates during freezing run
- Control multiple freezers with one instance of DeltaT
- Profiles with up to 32 different stages can be created
- Remote access of freezing run
- Password protected and tamper proof files
- Helps compliance with 21 CFR Part 11

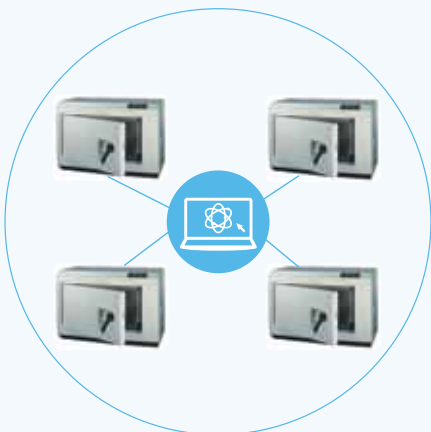



Freezer control with batch option DeltaT v7

	DeltaT v7
Create profiles	yes
Run Profiles	yes
View kryofiles	yes
Print kryofiles	yes
Export kryofiles	yes
Calibration	yes
Run qualification	yes
Batch mode	yes
Network capabilities	yes
Compatibility with the MR6 controller	yes



Specifications may change without notice, third party trademarks acknowledged. Ai030/V1